

Provicol® QM Aesthetic

TRANSLUCENT, TEMPORARY LUTING CEMENT WITH CALCIUM HYDROXIDE, EUGENOL-FREE



Provicol® QM Aesthetic

TRANSLUCENT, TEMPORARY LUTING CEMENT WITH CALCIUM HYDROXIDE, EUGENOL-FREE


Indications

Temporary luting of provisional or definitive inlays, onlays, crowns, partial crowns / veneers and bridges

Temporary filling of small one-surface cavities

Advantages

- Highly aesthetic, translucent appearance
- Increased strength for more retention and secure adhesion
- Optimum flowability for easy cementation
- Pinpoint stability allows precise removal of excess material
- · Accurate fit due to very low film thickness
- Eugenol-free: no inhibition of curing reaction of resinbased filling or luting materials
- Free from colophony: low allergy potential
- · Contains calcium hydroxide and zinc oxide
 - Calcium hydroxide supports the formation of tertiary dentine und has bacteriostatic effect
 - Zinc oxide is known for its antibacterial effect
- Radiopaque

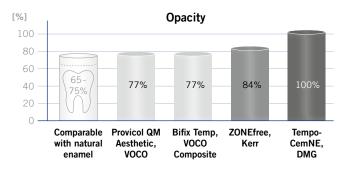

Provicol QM Aesthetic reaches a translucency, which is usually only reached by resin-based products. At the same time, Provicol QM Aesthetic combines all advantages of a luting cement and is compatible with a wide range of temporary and permanent restauration materials.

Source: internal measurement

Comparison translucency



Visual comparison of translucency and shade – Provicol QM Aesthetic also meets the demands in the anterior region, which were previously reserved for resin-based products.

Presentations

REF 1108 QuickMix syringe 5 ml, accessories REF 1109 Set II QuickMix syringe 3×5 ml $(1 \times \text{Provicol QM, } 1 \times \text{Provicol QM Plus,} 1 \times \text{Provicol QM Aesthetic}), accessories$

REF 2189 Mixing tips type 10, 50 pcs.

