


# SIMPLICITY IN 1 SHADE

## + X-TRA

- x-tra low shrinkage
- x-tra depth of cure
- x-tra high biocompatibility

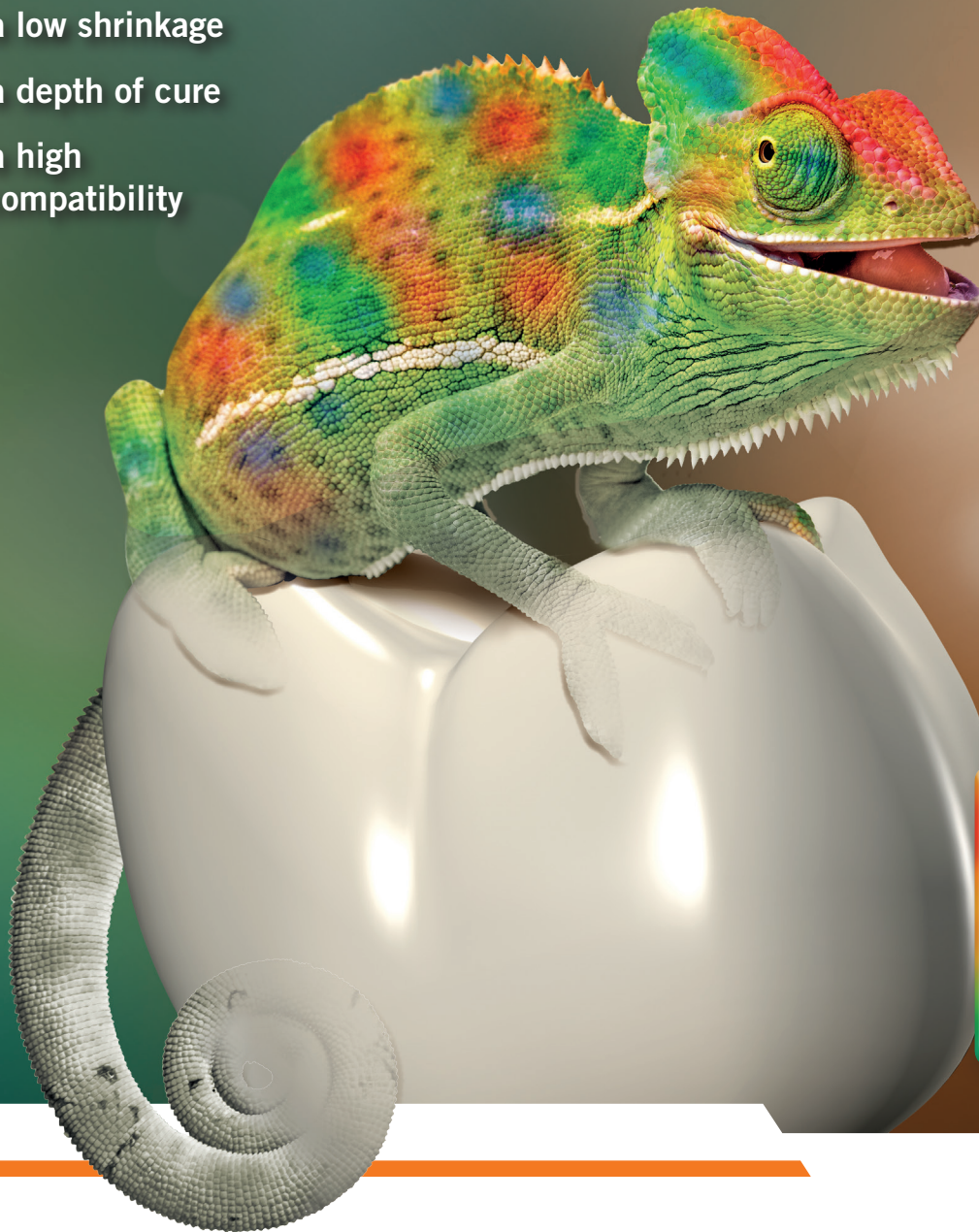

Visit  
[vocoamerica.com](http://vocoamerica.com)  
for your  
**FREE**  
Sample

## Admira Fusion x-tra

SINGLE-SHADE OMNI-CHROMATIC  
NANO-ORMOCER RESTORATIVE



# Admira Fusion x-tra

## SINGLE-SHADE OMNI-CHROMATIC NANO-ORMOCER RESTORATIVE

### Achieving exceptional shade matching through unique design, composition and formulation

The “nano” aspect of nano-ORMOCER technology refers to the Nano-hybrid filler technology utilized by Admira Fusion x-tra. These nano-hybrid fillers are composed of 60 % nano-particulate and 40 % micro and macro particulate – or glass filler immersed in an ORMOCER matrix. When the nano-particulate is viewed under an electron microscope, it is spherical in shape and measures 20–40 nanometers in diameter. It is this size and shape that dictates its interaction with light as it passes through the material. Due to its size and shape, Admira Fusion x-tra does not diffract or refract light. Rather, light passes through the nano-particulate, hits the surrounding tooth structure and returns to the human eye now influenced by the shade and color of that surrounding tooth structure. It is for this reason that Admira Fusion x-tra has the unique ability to match the entire shade range of human dentition in a single omni-chromatic universal shade.

### Advantages


- Only one omni-chromatic shade: no guess work and no more wasted shades
- 4 mm depth of cure – saves time
- Easy to polish, highly stain resistant with a high gloss
- Exceptional longevity with low initial shrinkage of only 1.25 % (v/v)
- 100 % BPA Free and no classic monomers for higher biocompatibility

### Indications

- Class I and II posterior restorations
- Base in Class I and II cavities
- Class V restorations
- Locking, splinting of loose anteriors
- Repairing veneers, small enamel defects and temporary C&B-materials
- Restoration of deciduous teeth
- Core build-up

### Presentation

- REF 2810 Syringe 3 g universal
- REF 2811 Caps 15 x 0.2 g universal

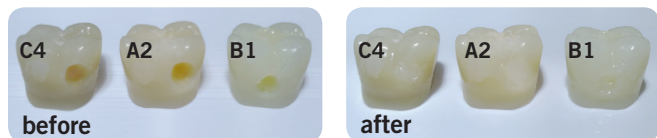

VOCO America Inc.  
1245 Rosemont Drive  
Suite 140  
Indian Land · SC 29707

www.vocoamerica.com  
infousa@voco.com

Toll-free phone: 1-888-658-2584  
Fax: 1-888-849-3989

